

# HDF5 MCS Raw Data Definition

Aus Multi Channel Systems Wiki

## Inhaltsverzeichnis

- 1 Definition of the HDF5 format for raw data
  - 1.1 Changelog
- 2 Hierarchy
- 3 Root-Folder "/"
- 4 Folder "Data"
- 5 Folder "Recording\_x"
- 6 Folder "AnalogStream"
- 7 Sub-folder "Stream\_x" of "AnalogStream"
- 8 Folder "FrameStream"
- 9 Subfolder "Stream\_x" of "FrameStream"
- 10 Subfolder "FrameDataEntity\_x"
- 11 Folder "EventStream"
- 12 Subfolder "Stream\_x" of "EventStream"
- 13 Folder "SegmentStream"
- 14 Subfolder "Stream\_x" of "SegmentStream"
- 15 **DataSubType-Average**: Subfolder "Stream\_x" of "SegmentStream"
- 16 Folder "TimeStampStream"
- 17 Subfolder "Stream\_x" of "TimeStampStream"
- 18 Comment

## Definition of the HDF5 format for raw data

MCS-HDF5 Protocol Type: **RawData** (Raw-Data protocol)

Protocol Version: **3** based on the definitions of RawDataFileIO in version 10.

All strings are only ASCII-encoded

## Changelog

Version 1:

- Initial draft


Version 2:

- New Root-Folder attributes added to detect name and version of the creating application and library

Version 3:

- Data structures for **DataSubType::Average** of **StreamType::Segment** added

## Hierarchy

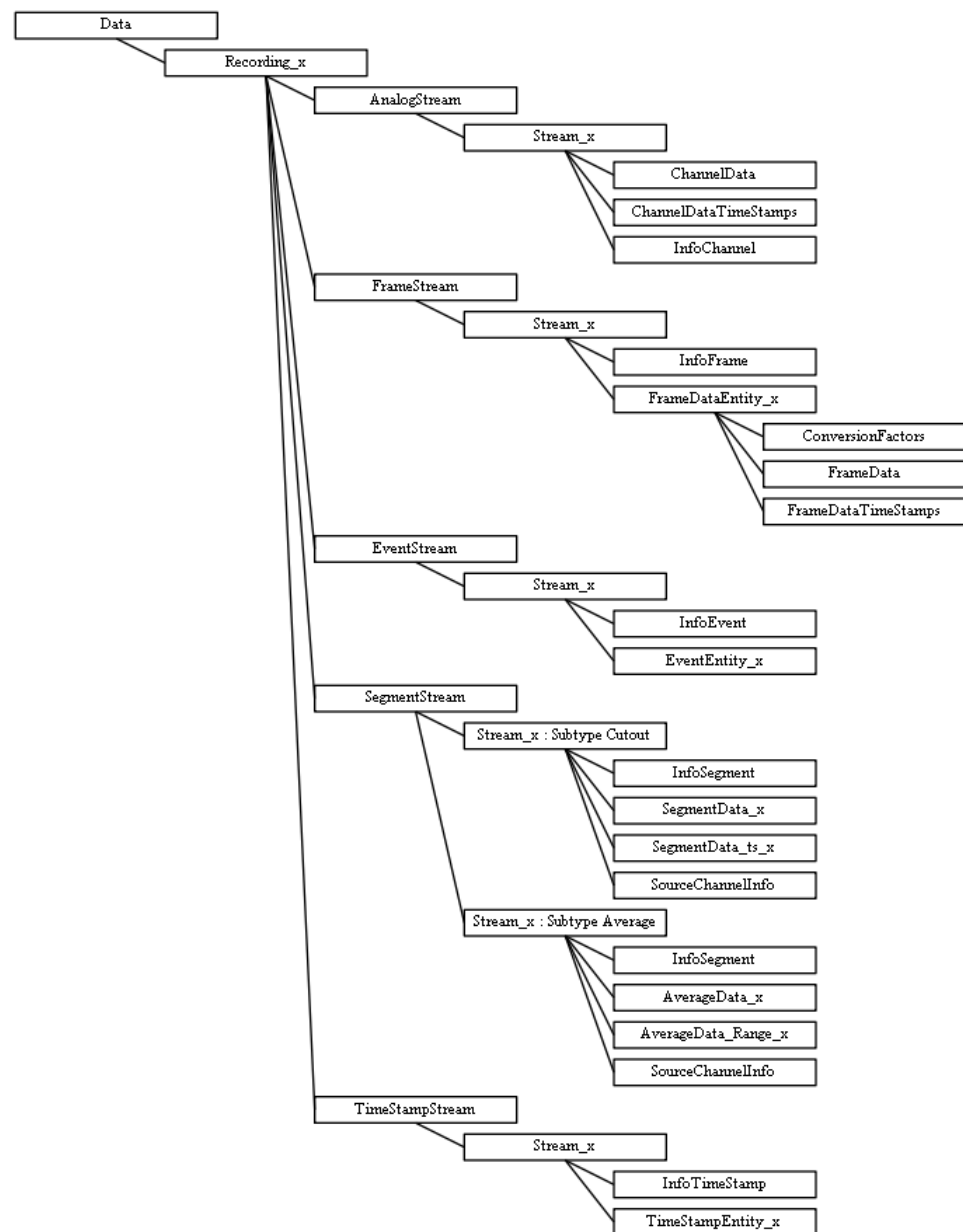

## Root-Folder "/"

Contains all information for one experiment - measured data (inside the folder **Data**) and a description (possibly in the future) inside the folder **Experiment/Description/...**

### Attributes:

Name	Description	Data Type	MCS-HDF5 Protocol Version
<a href="#">McsHdf5ProtocolType</a>	Type of the used MCS-HDF5 protocol definition (e.g. <b>RawData</b> for the raw data MCS-HDF5 definitions)	[String,Scalar]	1 ≤
<a href="#">McsHdf5ProtocolVersion</a>	Version number of the used MCS-HDF5 protocol	[Integer,Scalar]	1 ≤
<a href="#">GeneratingApplicationName</a>	Name of the application that generated this HDF5 file	[String,Scalar]	2 ≤
<a href="#">GeneratingApplicationVersion</a>	Version of the application that generated this HDF5 file	[String,Scalar]	2 ≤
<a href="#">McsDataToolsVersion</a>	Version of the McsDataTools library that was used by the application to create the HDF5 file	[String,Scalar]	2 ≤

**Datasets:**

- none

**Folder "Data"**

Navigation: /Data

Contains all recordings for this experiment.

**Attributes:**

Name	Description	Data Type
<a href="#">ProgramName</a>	Name of the recording program	[String,Scalar]
<a href="#">ProgramVersion</a>	Version number of the recording program	[String,Scalar]
<a href="#">MeaName</a>	Name of the recorded MEA	[String,Scalar]
<a href="#">MeaLayout</a>	Layout descriptor	[String,Scalar]
<a href="#">MeaSN</a>	Serial number of the MEA	[String,Scalar]
<a href="#">Date</a>	Date of the recording	[String,Scalar]
<a href="#">DateInTicks</a>	Date of the recording in .NET ticks (100 ns)	[Long(64-bit Integer),Scalar]
<a href="#">FileGUID</a>	GUID of the converted raw data file	[String,Scalar]
<a href="#">Comment</a>	Comment	[String,Scalar]

**Datasets:**

- none

**Folder "Recording\_x"**

Navigation: /Data/Recording\_x

Contains all recorded streams for recording x.

**Attributes:**

Name	Description	Data Type
<a href="#">RecordingID</a>	Recording ID	[Integer(32-bit Integer),Scalar]

<a href="#">RecordingType</a>	Recording type	[ <a href="#">String</a> , <a href="#">Scalar</a> ]
<a href="#">TimeStamp</a>	Start time of the recording in microseconds	[ <a href="#">Long</a> (64-bit Integer), <a href="#">Scalar</a> ]
<a href="#">Duration</a>	Total recording duration in microseconds (This duration can differ from the actual duration of the recorded data!!!)	[ <a href="#">Long</a> (64-bit Integer), <a href="#">Scalar</a> ]
<a href="#">Label</a>	Label	[ <a href="#">String</a> , <a href="#">Scalar</a> ]
<a href="#">Comment</a>	Comment	[ <a href="#">String</a> , <a href="#">Scalar</a> ]

**Datasets:**

- none

**Folder "AnalogStream"**

Navigation: /Data/Recording\_x/AnalogStream

(Organisational) folder for all channel-based streams of this recording

**Attributes:**

- none

**Datasets:**

- none

**Sub-folder "Stream\_x" of "AnalogStream"**

Navigation: /Data/Recording\_x/AnalogStream/Stream\_x

Container for an analog stream

**Attributes:**

Name	Description	Data Type	StreamInfoVersion
<a href="#">StreamInfoVersion</a>	Version number of the meta information structure	[ <a href="#">Int</a> (32-bit Integer), <a href="#">Scalar</a> ]	1 ≤
<a href="#">Label</a>	Label	[ <a href="#">String</a> , <a href="#">Scalar</a> ]	1 ≤
<a href="#">SourceStreamGUID</a>	GUID of the source streams	[ <a href="#">String</a> , <a href="#">Scalar</a> ]	1 ≤
<a href="#">StreamGUID</a>	GUID	[ <a href="#">String</a> , <a href="#">Scalar</a> ]	1 ≤
<a href="#">StreamType</a>	Type of the stream, e.g. <b>Electrode</b>	[ <a href="#">String</a> , <a href="#">Scalar</a> ]	1 ≤
<a href="#">DataSubType</a>	Sub-type of the analog stream (e.g. <b>Analog</b> )	[ <a href="#">String</a> , <a href="#">Scalar</a> ]	1 ≤

**Datasets:**

- Matrix [InfoChannel](#) →  $n \times 16$  matrix of describing information vectors for the n channels:
  - **Attributes:** [InfoVersion](#) → Version number of the Info-Objects [[Int](#)(32-bit Integer),[Scalar](#)]

Name	Description	Data Type	InfoVersion
<a href="#">ChannelID</a>	ID of the channel as given by the recording software	[ <a href="#">Int</a> (32-bit Integer), <a href="#">Array</a> (Size 1)]	1 ≤

RowIndex	Row number of this channel inside the <b>ChannelData</b> matrix where the data of this channel is stored	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$
GroupID	ID of the group that this channel belongs to	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$
Label	Label of the channel	[String,Array]	$1 \leq$
RawDataType	Type of the raw data	[String,Array]	$1 \leq$
Unit	Physical unit of the measured sensor value	[String,Array]	$1 \leq$
Exponent	Exponent $n \Rightarrow 10^n$ in which the channel values magnitude is measured (e.g. k,m, $\mu$ ,...)	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$
ADZero	ADC-Step that represents the 0-point of the measuring range of the ADC	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$
Tick	Sample tick $\Delta$ between two sample points of a channel in $\mu s$ $\Rightarrow$ sampling frequency = $1000000 / \Delta$	[Long(64-bit Integer),Array(Size 1)]	$1 \leq$
ConversionFactor	Conversion factor for the mapping ADC-Step $\Rightarrow$ measured value	[Long(64-bit Integer),Array(Size 1)]	$1 \leq$
ADCBits	Number of bits used by the AD-Converter	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$
HighPassFilterType	Type of the high-pass filter (empty string if not available)	[String,Scalar]	$1 \leq$
HighPassFilterCutOffFrequency	Cut-off frequency of the high-pass filter ('-1'-String if not available)	[String,Scalar]	$1 \leq$
HighPassFilterOrder	Order of the high-pass filter (-1 if not available)	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$
LowPassFilterType	Type of the low-pass filter (empty string if not available)	[String,Scalar]	$1 \leq$
LowPassFilterCutOffFrequency	Cut-off frequency of the low-pass filter ('-1'-String if not available)	[String,Scalar]	$1 \leq$
LowPassFilterOrder	Order of the low-pass filter (-1 if not available)	[Int(32-bit Integer),Array(Size 1)]	$1 \leq$

- 2-dimensional Data-Matrix **ChannelData**  $\rightarrow$  Data for sampled channels organized as  $n \times m$  matrix  $\Rightarrow$  one row per channel and one column per sample time point
  - reconstruct the value of the measured signal:  

$$y(\text{channel}, t_{ind}) = (\text{ChannelData}[\text{InfoChannel}[\text{channel}]. \text{RowIndex}, t_{ind}] - \text{ADZero}) * \text{InfoChannel}[\text{channel}]. \text{ConversionFactor} * 10^{\text{InfoChannel}[\text{channel}]. \text{Exponent}} \text{ in } \text{InfoChannel}[\text{channel}]. \text{Unit}$$
  - reconstruct the sample time point:  $t = t_{ind} * \text{InfoChannel}[\text{channel}]. \text{Tick}$  in  $\mu s$
- Matrix **ChannelDataTimeStamps**  $\rightarrow k \times 3$  matrix of segments where the rows are one segment and the columns are:
  - first column  $\rightarrow$  time stamp of the first sample point of the segment
  - second column  $\rightarrow$  first index (column) of the segment in **ChannelData**
  - third column  $\rightarrow$  last index (column) of the segment in **ChannelData**

## Folder "FrameStream"

Navigation: /Data/Recording\_x/FrameStream

(Organisational) folder for all frame-based streams of this recording

### Attributes:

- none

### Datasets:

- none

## Subfolder "Stream\_x" of "FrameStream"

Navigation: /Data/Recording\_x/FrameStream/Stream\_x

Folder that contains all Frame-Entities of one Frame-Stream:

#### Attributes:

Name	Description	Data Type	StreamInfoVersion
<a href="#">StreamInfoVersion</a>	Version number of the meta information structure	[ <b>Int</b> (32-bit Integer), <b>Scalar</b> ]	1 ≤
<a href="#">Label</a>	Label	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">SourceStreamGUID</a>	GUID of the source stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamGUID</a>	GUID	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamType</a>	Type of the stream <b>Frame</b>	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">DataSubType</a>	Sub-type of the event stream (e.g. <b>SpikeTimeStamp</b> )	[ <b>String</b> , <b>Scalar</b> ]	1 ≤

#### Datasets:

- Matrix [InfoFrame](#) →  $n \times 24$  matrix of describing information vectors for the  $n$  Frame-Entities:
  - Attributes:** [InfoVersion](#) → Version number of the Info-Objects [**Int**(32-bit Integer),**Scalar**]

Name	Description	Data Type	InfoVersion
<a href="#">FrameID</a>	ID of the frame entity as given by the recording software	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">FrameDataID</a>	ID of the frame entity inside the stream folder that maps this information vector to the entity folder ( <b>FrameDataID</b> → subfolder <b>FrameDataEntity_FrameDataID</b> )	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">GroupID</a>	ID of the group that this frame entity belongs to	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">Label</a>	Label of the entity	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">RawDataType</a>	Type of the raw data	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">Unit</a>	Physical unit of the measured sensor value	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">Exponent</a>	Exponent $n \Rightarrow 10^n$ in which the sensor values magnitude is measured (e.g. k,m, $\mu$ ,...)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">ADZero</a>	ADC-Step that represents the 0-point of the measuring range of the ADC	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">ADCBits</a>	Number of bits used by the AD-Converter	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">Tick</a>	Sample tick $\Delta$ between two frames in $\mu$ s $\Rightarrow$ sampling frequency = $1000000 / \Delta$	[ <b>Long</b> (64-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">HighPassFilterType</a>	Type of the high-pass filter (empty string if not available)	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">HighPassFilterCutOffFrequency</a>	Cut-off frequency of the high-pass filter ('-1'-String if not available)	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">HighPassFilterOrder</a>	Order of the high-pass filter (-1 if not available)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">LowPassFilterType</a>	Type of the low-pass filter (empty string if not available)	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">LowPassFilterCutOffFrequency</a>	Cut-off frequency of the low-pass filter ('-1'-String if not available)	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">LowPassFilterOrder</a>	Order of the low-pass filter (-1 if not available)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">SensorSpacing</a>	Distance between adjacent sensors in $\mu$ m	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">FrameLeft</a>	Sensor count of the left edge of the entity frame based on the reference frame	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤

<a href="#">FrameTop</a>	Sensor count of the top edge of the entity frame based on the reference frame	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$
<a href="#">FrameRight</a>	Sensor count of the right edge of the entity frame based on the reference frame	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$
<a href="#">FrameBottom</a>	Sensor count of the bottom edge of the entity frame based on the reference frame	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$
<a href="#">ReferenceFrameLeft</a>	Sensor count of the left edge of the reference frame (defined by the used sensor array)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$
<a href="#">ReferenceFrameTop</a>	Sensor count of the left edge of the reference frame (defined by the used sensor array)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$
<a href="#">ReferenceFrameRight</a>	Sensor count of the left edge of the reference frame (defined by the used sensor array)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$
<a href="#">ReferenceFrameBottom</a>	Sensor count of the left edge of the reference frame (defined by the used sensor array)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	$1 \leq$

## Subfolder "FrameDataEntity\_x"

Navigation: /Data/Recording\_x/FrameStream/Stream\_x/FrameDataEntity\_x

Contains all datasets of the Frame-Entity x

### Datasets:

- Matrix [ConversionFactors](#)  $\rightarrow n \times m$  matrix of conversion factors for the sensor array
- 3-dimensional Data-Cube [FrameData](#)  $\rightarrow$  cube of the frame data organized as one frame to one sample time point ( $n \times m$  matrix of sampled signal values per sensor)  $\times$  sample time points
  - reconstruct the value of the measured signal:  $y = (\text{FrameData}[x,y,t] - \text{ADZero}) * \text{ConversionFactors}[x,y]$
  - reconstruct the sample time point:
- Matrix [FrameDataTimeStamps](#)  $\rightarrow k \times 3$  matrix of segments where the rows are one segment and the columns are:
  - first column  $\rightarrow$  time stamp of the first sample point of the segment
  - second column  $\rightarrow$  first index (z-axis) of the segment in **FrameData**
  - third column  $\rightarrow$  last index (z-axis) of the segment in **FrameData**

### Datasets:

- none

## Folder "EventStream"

Navigation: /Data/Recording\_x/EventStream

(Organisational) folder for all event-based streams of this recording

### Attributes:

- none

### Datasets:

- none

## Subfolder "Stream\_x" of "EventStream"

Navigation: /Data/Recording\_x/EventStream/Stream\_x

Folder that contains all Event-Entities of one Event-Stream:

#### Attributes:

Name	Description	Data Type	StreamInfoVersion
<a href="#">StreamInfoVersion</a>	Version number of the meta information structure	[ <b>Int</b> (32-bit Integer), <b>Scalar</b> ]	1 ≤
<a href="#">Label</a>	Label	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">SourceStreamGUID</a>	GUID of the source stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamGUID</a>	GUID of the current stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamType</a>	Type of the stream <b>Event</b>	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">DataSubType</a>	Sub-type of the event stream (e.g. <b>StgSideband</b> , <b>UserInput</b> , <b>DigitalPort</b> )	[ <b>String</b> , <b>Scalar</b> ]	1 ≤

Sub-type Description:

- **StgSideband** → The event is associated to a STG sideband change.
- **UserInput** → The event is associated with an user input.
- **DigitalPort** → The event is associated with a digital port change.

#### Datasets:

- Matrix [InfoEvent](#) →  $n \times 7$  matrix of describing information vectors for the  $n$  Event-Entities:
  - **Attributes:** [InfoVersion](#) → Version number of the Info-Objects [**Int**(32-bit Integer),**Scalar**]

Name	Description	Data Type	InfoVersion
<a href="#">EventID</a>	ID of the event entity	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">GroupID</a>	ID of the group that the entity belongs to	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">Label</a>	Label of the entity	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">RawDataType</a>	Type of the raw data	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">RawDataBytes</a>	Number of bytes of the raw data type	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">SourceChannelIDs</a>	Comma separated list of ID's of (source) channel that were involved in the generation of this event	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">SourceChannelLabels</a>	Comma separated list of labels of the source channels	[ <b>String</b> , <b>Scalar</b> ]	1 ≤

- 2-dimensional matrix [EventEntity\\_x](#) →  $5 \times n$  matrix ⇒  $n$  events with describing vector (time stamp of event, duration of event, event info type, info 1, info 2)
  - **Attributes:** Short description of content
  - $t_{\text{event } i} = \text{EventEntity\_x}[0, i]$  in  $\mu s$
  - $\Delta_{\text{event } i} = \text{EventEntity\_x}[1, i]$  in  $\mu s$

## Folder "SegmentStream"

Navigation: /Data/Recording\_x/SegmentStream

(Organisational) folder for all segment-based streams of this recording. A segment is a cutout of parts of the sampled signal relative to an event, defined by a pre- and post interval.

#### Attributes:

- none



**Datasets:**

- none

**Subfolder "Stream\_x" of "SegmentStream"**

Navigation: /Data/Recording\_x/SegmentStream/Stream\_x

Folder that contains all Segment-Entities of one Segment-Stream:

**Attributes:**

Name	Description	Data Type	StreamInfoVersion
<a href="#">StreamInfoVersion</a>	Version number of the meta information structure	[ <b>Int</b> (32-bit Integer), <b>Scalar</b> ]	1 ≤
<a href="#">Label</a>	Label	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">SourceStreamGUID</a>	GUID of the source stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamGUID</a>	GUID of the current stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamType</a>	Type of the stream <b>Segment</b>	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">DataSubType</a>	Sub-type of the segment stream (e.g. <b>Spike</b> )	[ <b>String</b> , <b>Scalar</b> ]	1 ≤

**Datasets:**

- Matrix [InfoSegment](#) → n × 7 matrix of describing information vectors for the n Segment-Entities:
  - **Attributes:** [InfoVersion](#) → Version number of the Info-Objects [**Int**(32-bit Integer),**Scalar**]

Name	Description	Data Type	InfoVersion
<a href="#">SegmentID</a>	ID of the segment entity	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">GroupID</a>	ID of the group that the segment entity belongs to	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">Label</a>	Label of the entity	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">PreInterval</a>	Time interval in μs before the segment defining event occurred - definition of the beginning of the segment	[ <b>Int</b> (64-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">PostInterval</a>	Time interval in μs after the segment defining event occurred - definition of the end of the segment length of the segment = PreInterval + PostInterval in μs	[ <b>Int</b> (64-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">SegmentType</a>	Type of the segment (e.g. <b>SpikeCutout</b> )	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">SourceChannelIDs</a>	Comma separated list of ID's of (source) channels that the segments are taken from → Link to the <a href="#">SourceChannelInfo</a> matrix	[ <b>String</b> , <b>Array</b> ]	1 ≤

- 2-dimensional matrix [SourceChannelInfo](#) → n × 15 matrix ⇒ n of describing vectors for the n source channels, the structure is the same as in [ChannelInfo](#) used in section Sub-folder "Stream\_x" of "AnalogStream"
  - **Attributes:** [InfoVersion](#) → Version number of the Info-Objects [**Int**(32-bit Integer),**Scalar**]
- Vector [SegmentData\\_ts\\_x](#) → n time stamps in μs of the event triggering the segment, one for each of the n segments contained by segment entity x
- 2-dimensional matrix or 3-dimensional cube [SegmentData\\_x](#) → k × n matrix (k sample points for one segment, n number of sampled segments) or k × m × n cube (k sample points for one segment, m number of segments for one time point/for one multi-segment, n number of sampled multi-segments) of segment data:
  - **Attributes:** [SourceChannelID](#) → Comma separated list of ID's of (source) channels that the segments are taken from [**String**,**Scalar**] (the same as in InfoSegment, repeated for clarification)
  - reconstruct the value of the measured segment signal (only one segment  $id_{\text{segment}}$  → 2-dimensional matrix **M**[row,col]):
    - $t_{ind}[row, col] = \text{SegmentData}_{ts\_x}[col] + (row - 1) * tick_{\text{source-channel}} - \text{PreInterval}$  in μs
    - $y(id_{\text{segment}}, t_{ind}(row, col)) = (\text{SegmentData}_{x[ow, col]} - \text{ADZero}_{\text{source-channel}}) * \text{ConversionFactor}_{\text{source-channel}} * 10^{\text{Exponent}_{\text{source-channel}}}$  in InfoChannel[source-channel]. Unit

- reconstruct the value of the measured segment signal (m segments  $\rightarrow$  multi-segments  $\rightarrow$  3-dimensional cube  $\mathbf{M}[\text{row}, \text{col}, \text{z}]$ ):
  - $\text{col} \rightarrow id_{\text{segment}} \rightarrow \text{source-channel}$
  - $t_{\text{ind}}[\text{row}, \text{col}, \text{z}] = \text{SegmentData ts x}[\text{z}] + (\text{row} - 1) * \text{tick}_{\text{source-channel}[\text{col}]} \text{ in } \mu\text{s}$
  - $y(id_{\text{segment}}, t_{\text{ind}}(\text{row}, \text{z})) = (\text{SegmentData x}[\text{row}, \text{col}, \text{z}] - \text{ADZero}_{\text{source-channel}[\text{col}]}) * \text{ConversionFactor}_{\text{source-channel}[\text{col}]} * 10^{\text{Exponent}_{\text{source-channel}[\text{col}]}} \text{ in } \text{InfoChannel}[\text{source-channel}[\text{col}]]. \text{Unit}$

## DataSubType-Average: Subfolder "Stream\_x" of "SegmentStream"

Navigation: /Data/Recording\_x/SegmentStream/Stream\_x

Folder that contains all Segment-Entities of one Segment-Stream with **DataSybType == Average**:

**Attributes:** no difference to the standard case above

**Datasets:**

- Matrix [InfoSegment](#): no difference to the standard case above
- Matrix [SourceChannelInfo](#): no difference to the standard case above
- (3  $\times$  n) matrix [AverageData\\_Range\\_x](#)  $\rightarrow$  (**start, end, count**) per segment average  $\times$  count of segment averages contained by segment entity x. **start** and **end** denote the start and end timestamp in  $\mu\text{s}$  of the interval that contains all averaged segments. **count** is the number of averaged segments.
  - Attributes:** description of the content
- (2  $\times$  k  $\times$  n) cube [AverageData\\_x](#)  $\rightarrow$  (mean and standard deviation)  $\times$  k sample points of the segment  $\times$  n number of segment averages
  - Attributes:** description of the content
  - reconstruct the value of the mean and standard deviation of the average segment (n average segments  $\rightarrow$  3-dimensional cube  $\mathbf{M}[\text{row}, \text{col}, \text{z}]$ ):
    - row: mean  $\rightarrow$  row = 0; StdDev  $\rightarrow$  row = 1
    - col:  $t_{\text{ind}}(\text{col}) = (\text{col} - 1) * \text{tick}_{\text{source-channel}} \rightarrow \text{time range } (0, \text{PreInterval}[\text{SegmentID}] + \text{PreInterval}[\text{SegmentID}]) \text{ in } \mu\text{s}$
    - z:  $z = id_{\text{average}}$  (number of average segment)
    - $\text{Mean}(id_{\text{average}}, t_{\text{ind}}(\text{col})) = (\text{AverageData x}[0, \text{col}, id_{\text{average}}] - \text{ADZero}_{\text{source-channel}}) * \text{ConversionFactor}_{\text{source-channel}} * 10^{\text{Exponent}_{\text{source-channel}}} \text{ in } \text{InfoChannel}_{\text{source-channel}}. \text{Unit}$
    - $\text{StdDev}(id_{\text{average}}, t_{\text{ind}}(\text{col})) = \text{AverageData x}[1, \text{col}, id_{\text{average}}] * \text{ConversionFactor}_{\text{source-channel}} * 10^{\text{Exponent}_{\text{source-channel}}} \text{ in } \text{InfoChannel}_{\text{source-channel}}. \text{Unit}$

## Folder "TimeStampStream"

Navigation: /Data/Recording\_x/TimeStampStream

(Organisational) folder for all TimeStamp-based streams of this recording

**Attributes:**

- none

**Datasets:**

- none

## Subfolder "Stream\_x" of "TimeStampStream"

Navigation: /Data/Recording\_x/TimeStampStream/Stream\_x

Folder that contains all TimeStamp-Entities of one TimeStamp-Stream:

**Attributes:**

Name	Description	Data Type	StreamInfoVersion
<a href="#">StreamInfoVersion</a>	Version number of the meta information structure	[ <b>Int</b> (32-bit Integer), <b>Scalar</b> ]	1 ≤
<a href="#">Label</a>	Label	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">SourceStreamGUID</a>	GUID of the source stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamGUID</a>	GUID of the current stream	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">StreamType</a>	Type of the stream <b>TimeStamp</b>	[ <b>String</b> , <b>Scalar</b> ]	1 ≤
<a href="#">DataSubType</a>	Sub-type of the TimeStamp stream (e.g. <b>NeuralSpike</b> )	[ <b>String</b> , <b>Scalar</b> ]	1 ≤

Sub-type Description:

- **NeuralSpike** → The entity contains time stamps of neural spikes

**Datasets:**

- Matrix [InfoTimeStamp](#) →  $n \times 7$  matrix of describing information vectors for the  $n$  Event-Entities:
  - **Attributes:** [InfoVersion](#) → Version number of the Info-Objects [**Int**(32-bit Integer),**Scalar**]

Name	Description	Data Type	InfoVersion
<a href="#">TimeStampEntityID</a>	ID of the event entity	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">GroupID</a>	ID of the group that the entity belongs to	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">Label</a>	Label of the entity	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">Unit</a>	Physical unit of the measured sensor value	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">Exponent</a>	Exponent $n \Rightarrow 10^n$ in which the channel values magnitude is measured (e.g. k,m,μ,...)	[ <b>Int</b> (32-bit Integer), <b>Array</b> (Size 1)]	1 ≤
<a href="#">SourceChannelIDs</a>	Comma separated list of ID's of (source) channel that were involved in the generation of this event	[ <b>String</b> , <b>Array</b> ]	1 ≤
<a href="#">SourceChannelLabels</a>	Comma separated list of labels of the source channels	[ <b>String</b> , <b>Scalar</b> ]	1 ≤

- Vector [TimeStampEntity\\_x](#) →  $n$  time stamps in μs

**Comment**

All time-related information except dates (100ns ticks) are given in microsecond ticks!!

▪